

topolino capo 2K

The camera system with 2K resolution for complex templates



The camera system topolino capo 2K combines highest resolution (2,5x higher than HD) with a huge magnification range from approx. 0,9x to 100x in reading mode*. Texts and pictures even larger than DIN A3 are being displayed contrast enhanced and razor sharp. These characteristics are unique and thus ideal for professional use.

The integrated HDMI- and USB 3.0 interfaces make sure that topolino capo 2K can be connected to a 2K monitor, or a PC or notebook.

The optionally available battery and XY table upgrade comfortability of topolino capo 2K. For transport, fold up the unit easily and space savingly.

Your benefits at a glance

- transportable 2K reading-/ distance viewing camera in highest resolution
- up to 2048 x 1280p*
- approx. 22x lossless magnification in FullHD, magnification in reading mode
- approx. 0,9x to 100x*
- compatible with Windows 10, 11, macOS und iPad
- razor sharp display of more than a DIN A3 page
- overview mode via fading cross display
- HDMI- and USB 3.0 interface
- compatible with SuperNova and Zoomtext, e.g. for display in split screen mode
- intuitive and flexible operation: via camera, remote control or PC keyboard

*depending on monitor size and connection (USB 3.0 or HDMI)

topolino capo 2K

order code:
TOCA2K11



60 Hz camera, progressive scan, system resolution of up to 2048 x 1280 pixels in distance viewing mode
autofocus (can be switched ON/OFF)



approx. 0,9x to 100x
in document mode

*depending on monitor size and connection (USB 3.0 or HDMI)



true colors

semi colors: **F F F F F**

individually selectable semi color modes



weight approx. 3,9 kg

dimensions: (W x H x D)

approx. 31 x 58,5 x 46 cm,

folded approx. 31 x 15,5 x 46 cm



compatible with Windows 10 and 11,
macOS, iPad (with USB-C port)



connections: HDMI port, USB 3.0 port



scope of supply: remote control, power supply, USB 3.0 cable



available options: mobile XY table, XY table with friction brake;
PC, notebook, monitor,
PC magnification software

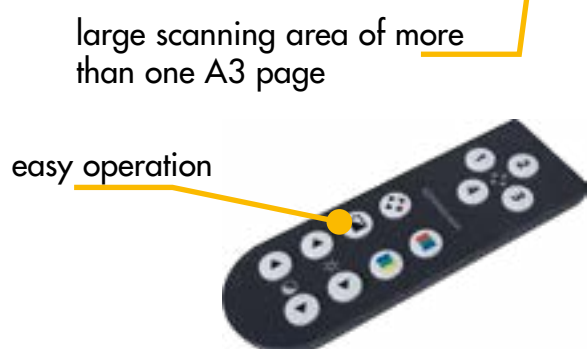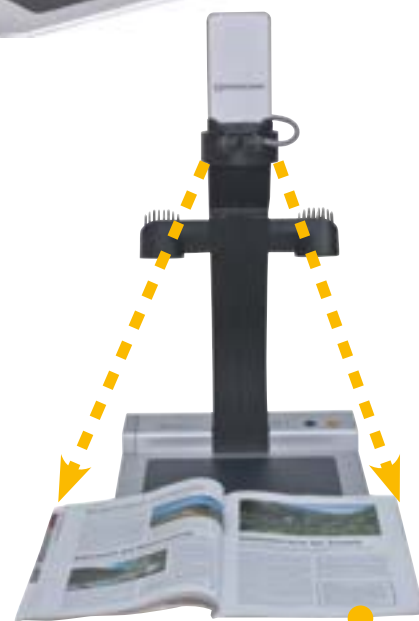

Reinecker Vision GmbH
Head office Germany
Sandwiesenstraße 19
D-64665 Alsbach-Hähnlein

Tel. +49 6257 9311-0

Fax +49 6257 9311-8008
international@reineckervision.de

Your authorized distributor: