


# Videomatic

The contemporary video magnifier

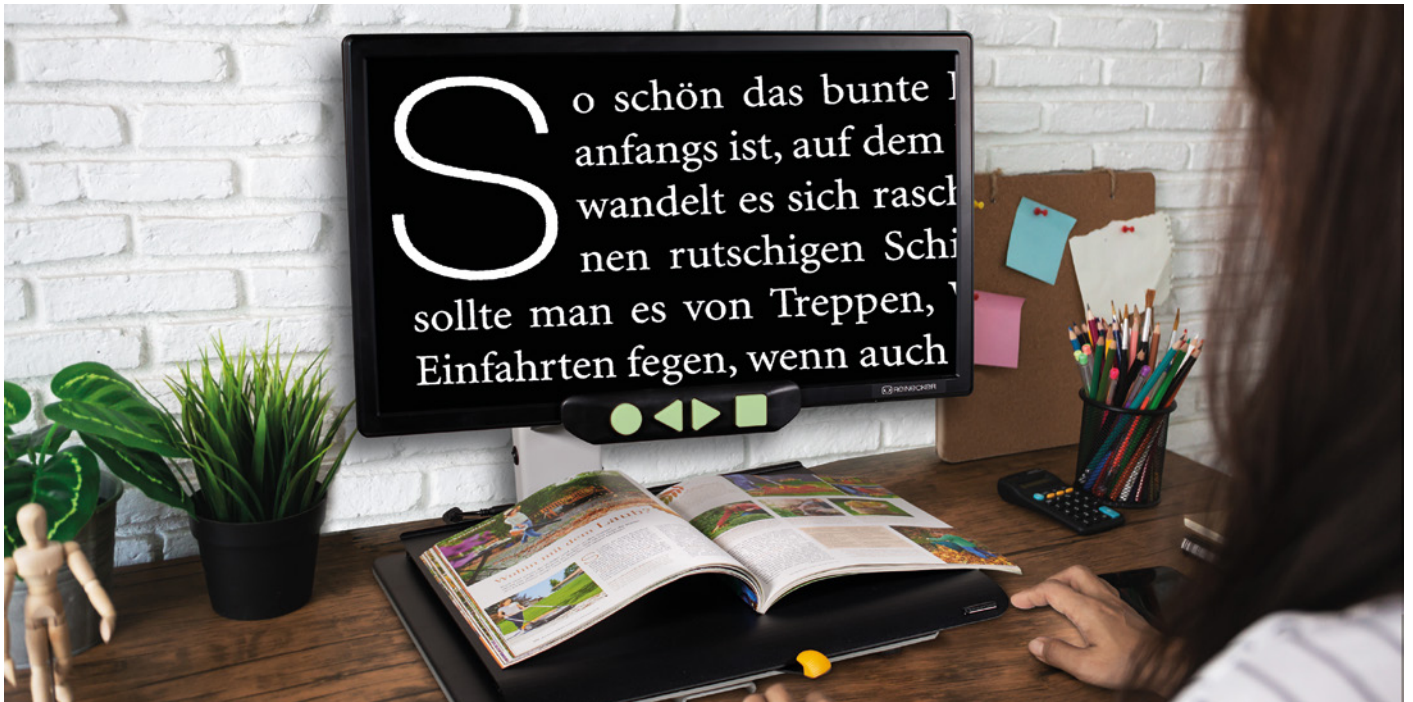

In 1971, the first video magnifier manufactured by Reinecker Vision, was named Videomatic. Ever since it stands for comfort, innovation and quality.

We are now presenting Videomatic in a new design and combining proven quality features with contemporary design and innovative technology.

The high-resolution flat screen is mounted on an ergonomic swivel-tilt arm that helps the user to position it ergonomically

according to individual needs. The operating panel is clearly structured with four large, tactile controls that make operation of Videomatic easy and intuitive.

Videomatic is equipped with a fullHD camera which magnifies texts of books or newspapers and other printed material such like letters. You can easily write or do cross words on the movable XY table with friction brake.

## Your benefits at a glance

- high resolution 20"-class LED flat screen
- ergonomic and comfortable positioning of the height adjustable and inclinable monitor
- three level light adjustment
- individual positioning of text or pictures underneath the camera due to movable XY table with friction brake
- ease of operation via four clearly structured and tactile controls
- overview mode: at the push of a button text will be zoomed out for better orientation at higher magnification levels; line masking and shadowing offer a targeted orientation on screen

## Videomatic



### LED flat screen

LED flat screen 20" class |  
approx. 49 cm  
adjustable in height, inclinable



### Focus camera:

FullHD camera with digital zoom



magnification approx. 4x to 42x \*

*\*refers to factory settings, further magnification levels to be adjusted within the internal menu*



true colors, semi colors: **F F F F F**  
individually selectable semi color modes




weight approx. 12 kg

dimensions approx. (W x H x D)  
46.5 x 54 x 49 cm



scope of supply: comfortable XY table,  
line masking, shadowing,  
overview mode,  
dimming function,  
individual menu  
available options: XY table

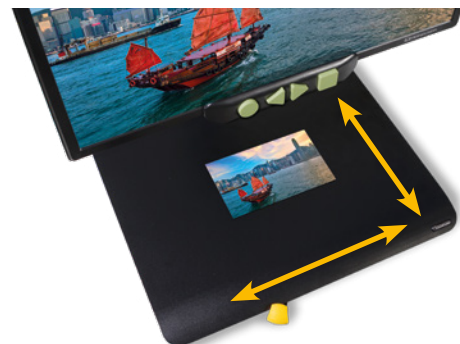

Reinecker Vision GmbH  
Head office Germany  
Sandwiesenstraße 19  
D-64665 Alsbach-Hähnlein

**Tel. +49 6257 9311-0**

Fax +49 6257 9311-8008  
international@reineckervision.de

Your authorized distributor: