

mezzo 2K

The battery-powered video magnifier with 2K document/distance camera for a high quality image



The combination of a high-resolution monitor and a hand-operated 2K camera with up to 22x lossless magnification ensures a top quality image even at high magnification levels.

The powerful rechargeable battery lasts up to 5 hours and guarantees a smooth workflow even if no power outlet is close by.

Thanks to its tiltable reading/distance camera, you can quickly switch between distance and reading mode and your last settings in terms of color contrast and magnification will always be saved. You can adjust the magnification either on the monitor's control panel or directly on the camera, e.g. while moving the distance camera.

Your benefits at a glance

- battery-powered and compact video magnifier with high resolution 16" monitor
- 2K reading/distance camera with highest resolution up to 2048x1152 pixels
- approx. 22x lossless magnification, magnification in document mode approx. 1.7 to 46x*
- dimmable lighting
- automatic saving of the last reading and distance mode settings
- movable monitor with spring compensation can be easily positioned and enables an ergonomic reading position
- can be set up and folded in a few simple steps, space-saving design

**refers to factory settings, further magnification levels to be adjusted within the internal menu*

mezzo 2K



LED flat screen

16" class | approx. 41 cm
(1920x1080 pixel)

Comfortable auto-focus
with one-touch deactivation



2K reading/distance camera:

up to 2048 x 1152 pixel



Magnification approx. 1.9x - 53x*
in document mode,

Lossless magnification up to
approx. 22x

**Refers to factory settings, further magnification
levels to be adjusted within the internal menu*



Battery operation:

Li-ion rechargeable battery

Operating time up to approx. 4.5 h
(depending on lighting)

Charging time approx. 5 h



Color representation:

true colors, semi colors:



individually selectable semi color modes



Weight approx. 5.9 kg

Dimensions approx. (W x H x D)

38 x 47 x 52 cm


Folded 38 x 19 x 52 cm



Scope of supply: carrying bag,
AC adapter, user manual



Available options: XY table

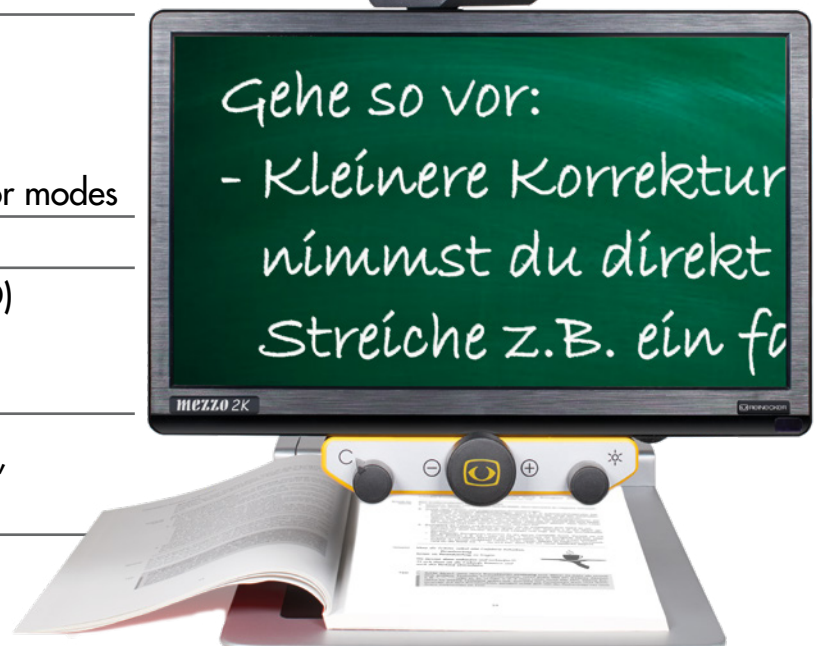

Reinecker Vision GmbH
Head office Germany
Sandwiesenstraße 19
D-64665 Alsbach-Hähnlein

Tel. +49 6257 9311-0

Fax +49 6257 9311-8008
international@reineckervision.de

Your authorized distributor: