

# MANO 6

The fullHD hand-held magnifier in wide-screen format



# We're here for you

We at Reinecker Vision are happy and proud to keep our traditional, independent, family owned standing since 1971 and always put our eye on intelligent and innovative product developments, all made in Germany.

With all of our product innovations that we export to more than 50 countries, we focus first on the demands of our low vision individuals.

Whether in the home, educational or professional environment we design products to help our customers reach their full potential.

Your benefits:

- almost 50 years of experience in low vision solutions
- short decision making progress
- individual support
- one face to the customer care

**Reinecker Vision**

# Everything from one source

If your vision declines considerably at an older age or as a result of an eye disease, a wide range of assistive devices can help you manage your everyday life. We are happy to help you put together individual solutions to improve your daily life.

## Our product range for your well-being

- ▶ **Electronic magnifiers: MANO**
- ▶ **Video magnifiers: mezzo, VEO, Videomatic**
- ▶ **Camera systems: topolino**
- ▶ **Reading systems: VOX**
- ▶ **Accessories: XY-tables, height adjustable tables, etc.**

## Product range



PRODUCT VIDEO



click to view



# MANO 6

MANO 6 is equipped with a high resolution touch display as well as with four tactile buttons for easy and comfortable operation. In addition, magnification level can be adjusted and the image can be scrolled by zoom gestures with your fingers. By means of the tactile buttons at the right and the left hand side of the display, different contrast modes, magnification levels can be adjusted and pictures can be taken.

With a height of only 9cm MANO 6 sits well in your hand. Due to the widescreen format, the reading flow is not interrupted because even at a high magnification level you can view several letters at the same time.

MANO 6 turns ON automatically when you fold out the stand. You can now place the magnifier on top of your reading template and start reading in an ergonomic reading angle.

Alternatively, you can also use the handle, e.g. to have a closer look at price tags while shopping. An upholstered bag protects MANO 6 from damages on the go.

**The fullHD hand-held magnifier in wide-screen format**



## MANO 6 – your benefits at a glance

- intuitive operation with only four buttons and touch display (diagonal of approx. 14cm)
- magnification adjustment either via tactile buttons or by zoom gestures on the touch display
- continuous magnification of approx. 2x – 35x
- equipped with an infinitely variable fold out handle which can be used as stand
- images can be saved in picture gallery and transferred to a micro SD-card
- artificial colors, line ruler, line masking, auto ON/OFF, button beep, etc.




**intuitive**




**comfortable**



**versatile**

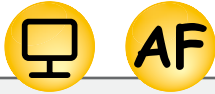

touch display with approx.  
14 cm diagonal  
snapshot function, memory of  
approx. 2GB; images can be  
transferred to a microSD card  
autofocus (can be switched off)

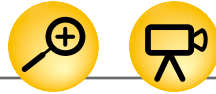

approx. 2x-35x  
system resolution  
1920 x 1080 pixel



Lithium-Polymer-battery,  
operating time approx. 3 h  
charging time approx. 2.5 h



scope of supply:  
carrying bag, short and long  
hand strap, power supply with  
country specific cord, micro fibre  
cloth, user manual, quick refer-  
ence guide



true colors,  
up to 16 selectable semi  
colors



weight approx. 274 g  
W x H x D approx.  
16 x 9,1 x 2,1 cm





# Reinecker USA LLC

145 River Rock Drive  
Buffalo, NY 14207  
USA

Tel. 716-348-3517

Fax 413-403-9777

Toll free: 866-733-2352

[sales@reineckerusa.com](mailto:sales@reineckerusa.com)

[www.reineckerusa.com](http://www.reineckerusa.com)